

BLACKHEAD REMOVER

Vacuum constriction · 3 Intensity · 3 replacement tips



CONTENTS



01		Product introduction
02		Operation
03		Blackhead remover Tips
04		Functions
05		Packaging
06		ODM&OEM
07		Certificatesions
08		Contact US

3D curved surface nozzle design

Smooth and delicate angle for better fit to the skin



3 Intensity

Low Intensity: Gentle adsorption for dry skin

Mid Intensity: Medium strength suction for mixed skin

High Intensity: Strong adsorption for oily skin

Original Design

Ergonomic design restores comfortable feeling,
Seamless design



Operation



Step 1

After cleaning, Steam face with hot towel or derived liquid for 3 to 5 minutes to open pores.



Step 2

Select suitable suction mode to suck out the white spots and dirt around the skin,
DO NOT stay on the sample place for 3 seconds



Step 3

Use the shrinking liquid to restore the treated pore. Tap gently until absorbed

Blackhead remover Tips



1.Oval Tip

For daily skin care and management of wrinkles at the corner of the mouth, nose or cheek. Use in a down to up motion for lifting facial skin.



2.Microcrystalline Tip

For effectively exfoliating dead skin for smooth, brighter looking skin and increasing skin surface absorption of your moisturizer, cream, and serum



3.Round Tip

For powerfully removing blackheads, acne, pimples, white comedone, grease, dirt and cosmetic residues. Recommended for use once a week, five minutes at a time.

Functions




Press the on/off button once: Low Intensity- Suitable for sensitive skin.
Press the on/off button twice: Mid Intensity- Suitable for areas where blackheads are not obvious.
Press the on/off button three times: High Intensity- Suitable for areas with large pores and obvious blackheads.

- The Red indicator light on when the battery power lower than 30%
- The Blue indicator light shows which intensity it's working
- The Green indicator light on means that it is fully charged



Packaging

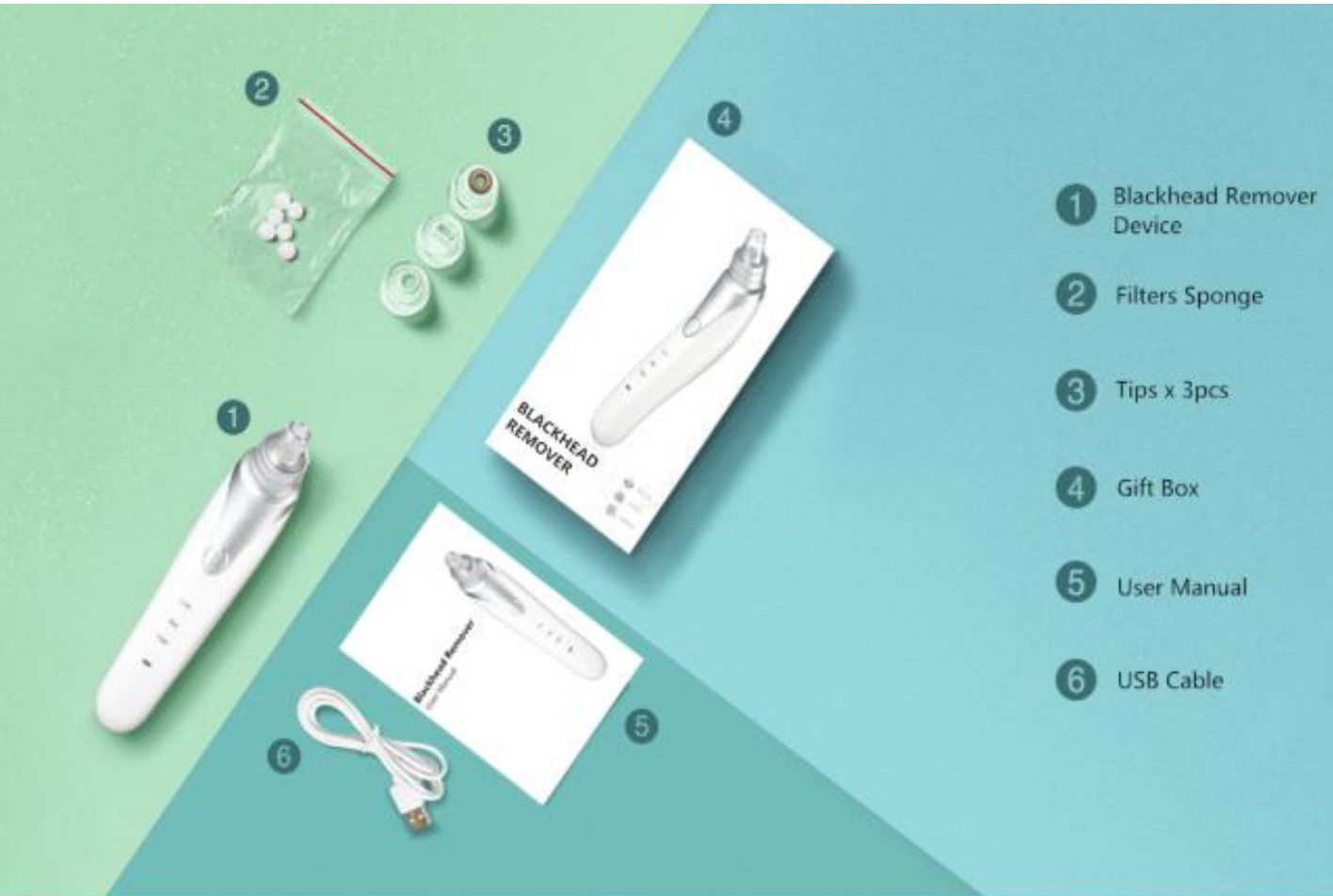

Unit pack Weight:	315g
Unit size:	195*41MM
Package Size:	105*206*43MM
Qty :	60pcs/ctn
Carton :	51*43*33cm
N.W./G.W.:	8.9/20.9KG

Customize Color



Customize Logo



Customize Box



Certifications



ROHS:BSTDG200313695601CC



FCC:BSTXD200313895602EC



CE:BSTXD200313695601EC

방송통신기자재등의 적합등록 필증 Registration of Broadcasting and Communication Equipments	
상호 또는 성명 Trade Name or Register	SHENZHEN BASELAB TECHNOLOGY LLC
기자재명칭(제품명칭) Equipment Name	BLACKHEAD REMOVER WITH HOT & COOL COMPRESS
기본모델명 Base Model Number	ALHTY07K-02
파생모델명 Series Model Number	ALHTY07K-06, ALHTY07K-04
등록번호 Registration No.	R-R-BTI-ALHTY07K-02
제조사/제조(조합)국가 Manufacturer/Country of Origin	SHENZHEN BASELAB TECHNOLOGY LLC / 중국
등록연월일 Date of Registration	2019-10-30
기타 Others	
<p>위 기자재는 「전자파」 제58조의2 제3항에 따라 등록되었음을 증명합니다. It is verified that foregoing equipment has been registered under the Clause 3, Article 58-2 of Radio Waves Act.</p> <p>2019년(Year) 10월(Month) 30일(Day)</p> <p>국립전파연구원 Director General of National Radio Research Agency</p> <p>※ 적합등록 방송통신기자재는 반드시 "적합성평가표시"를 부착하여 유통하여야 합니다. 의면시 카테르프 처분 및 등록이 취소될 수 있습니다.</p>	

KC:R-R-WKM-DSS

Contact US

SAY
Hello we'd love to talk...



QPP LIMITED

ADD: Building A, No. 6 Compound, Dongcai Industrial Zone, Gushu Area, Xixiang, Baoan District, Shenzhen, 518101, P.R.China

TEL: +86 021 6059 9416

Website: www.jannu.cn

Whatsapp/Viber/Wechat: +86 13129572930